• ИРРИГАТОР ПОЛОСТИ РТА • ORAL IRRIGATOR •
• MUNDDUSCHE • IRRIGADOR BUCAL • IRRIGATEUR BUCCAL •

МОДЕЛЬ/MODEL/MODELL/MODELO
MODÈLE: VIP-003

VES electric

RU РУКОВОДСТВО ПО ЭКСПЛУАТАЦИИ
EN INSTRUCTION MANUAL
DE BEDIENUNGSANLEITUNG
ES MANUAL DE INSTRUCCIONES
FR MODE D'EMPLOI

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EN DENTAL WATER JET

MODEL: VIP-003

Do not use this appliance for more than 5 minutes in each two-hour period. Longer usage can overheat the motor and cause appliance damage as well as void the warranty.

Dear customer! We highly appreciate your choosing VES Electric home appliances. From our side we guarantee excellent quality, perfect functioning and safety of your appliance upon keeping to all the rules of home exploitation.

Before using the appliance read the instruction manual carefully.


Technical data: rated power: 45 W, rated voltage: 220 - 240 V
Rated frequency: 50/60 Hz / Maximum pressure of the jet: 920 kPa

Recommendations for buyer


When buying the appliance, check it for the absence of mechanical damage of the construction, check the external design and packing, the functioning, the presence of all accessories in set. Note: operation checking - is the duty of the seller. Please check the warranty card and the accuracy of its filling (presence of the seller's seal). You must save the bill, instruction manual, warranty card until the end of the warranty period. In case of bill and warranty card loss they cannot be renewed. In the case of contacting the service center, if it is impossible to determine the date of product purchasing, the warranty period is determined from the manufacturing date.

DESCRIPTION

1. Stepless water pressure regulator
2. On/off switch / 3. Water tank
4. Water tank lid with tips storage compartment / 5. Jet tip
6. Rotary mechanism for the tip
7. Eject button for tip / 8. Grip
9. Pause button / 10. Water hose
11. Base

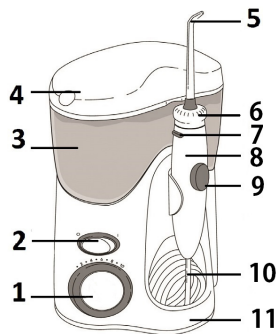ACCESSORIES

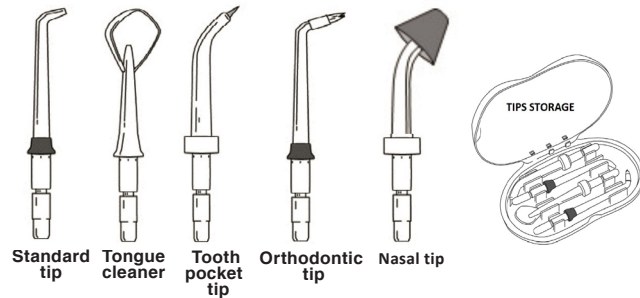BEFORE FIRST USE Remove the appliance from the package. Do not leave the packing (plastic bags, polystyrene etc.) within the reach of children since they could be dangerous. Keep the package for appliance storage.


MAINTENANCE Note: If your physician recommended you to take antibiotics before having oral hygiene procedures, consult your dentist before using given appliance or any other appliance for oral hygiene.

WATER JET PRETREATMENT Detach water tank with lid from the base. Fill water tank with warm water, close the lid. Place the water tank on the base, press it slightly.

Attention! Do not use bleaching or iodine-containing fluids.

Fill the water tank with water or fluids (rinsing or therapeutic fluids), recommended by your dentist and nothing else. The liquid should not leave sediment or be viscous, this may soil the appliance.

TIP INSTALLATION AND DETACHMENT Place the tip into the hole on the upper

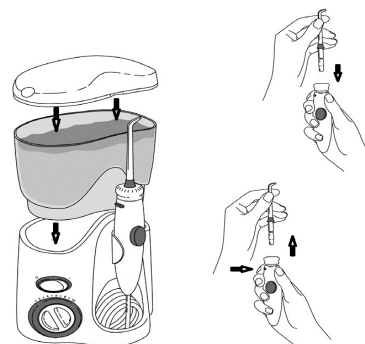surface of appliance grip. The color ring on the tip should closely touch the grip, in case of correct tip installation. To detach the tip, press tip eject button and pull the tip upwards.

WATER PRESSURE LEVEL SETTING Move water pressure regulator located on the base into the utmost left position – minimum water pressure. Slowly increase the pressure to the required level, you may choose it yourself or with the help of your dentist.

LOCATION IN THE MOUTH Locate the tip in the mouth. Direct the tip end at the teeth. Switch the appliance on. Set the required pressure level. Direct water jet, turning the tip. Close your lips to avoid water splashing, but do not close your mouth allowing water to freely flow away. For better result please start with side teeth and molars, moving it to frontal teeth. Move tip from tooth to tooth making pauses between moves. Go on cleaning until teeth surrounding zones and zones between teeth are clean.



OPERATION TIPS Caution! Do not direct water jet under your tongue, into the ear, nose and any body parts. Water jet pressure can cause serious harm to these body parts. Do not use this appliance if you have jewelries in the mouth; take them away before using the appliance. Do not use the appliance if you have open sore on tongue/ in the mouth.

PAUSE You can permanently stop water supply by pressing pause button on the appliance grip.


TIPS

JET TIP is used for regular teeth cleaning.

TONGUE CLEANER is used to remove bacteria and fur. Locate the tip in the mouth center carefully. Turn the appliance on, move water pressure regulator into the utmost left position and slightly pressing the tip pull it along the tongue in the lips direction. Repeat this process several times to clean the entire surface of the tongue. You can increase water pressure during the process if it feels comfortable.

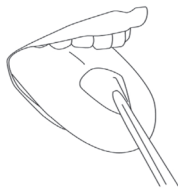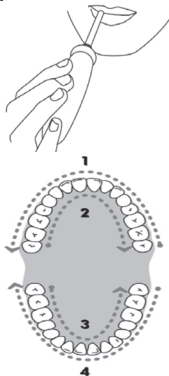
TOOTH POCKET TIP (periodontal) allows supplying water or antibacterial solvent directly into the teeth pockets. For this direct soft end of the tip at the tooth at 45 degrees angle and gently insert the tip under the gum line, into the tooth pocket. Turn the appliance on, move water pressure regulator into the utmost left position and perform cleaning moving the tip between the teeth and along the gum line.

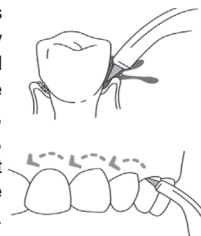
ORTHODONTIC TIP

Orthodontic tip is used for cleaning inaccessible places: around dental brackets, crown works, frameworks and dentures. For this direct the tip at the teeth at 90 degrees angle, making slight effort move the tip from molars to the frontal teeth clean the space along gum line, hold the tip between teeth for some time to clean the space between them and around dental bracket (crownwork), then clean inner side of the tooth and get to the next tooth.

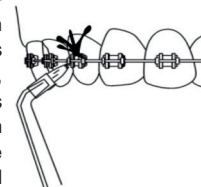NASAL TIP

Nasal tip is used for nose rinsing and ENT diseases treatment. Rinsing nose with salty water or water with medications, unblocks the nose, releases breathe, removes mucus from nasal sinuses. Use it in accordance with the below pictures.

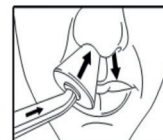CLEANING Unplug the appliance from mains. Detach water tank. Pour water remnants out of water tank. Rinse water tank with clean water, let dry. Wipe the base with wet cloth. Install water tank back onto the base. Roll power cord carefully. For appliance cleaning use only the clean water or neutral cleansers. Do not use the caustic or abrasive cleansers (e.g. vinegar or decalcifying agents) as they

may damage the appliance. Do not wash the tips with water hotter than 40 C.

- If you used the mouth wash. its remnants must be cleaned away. For this clean the tank with clean water, pour the filtered clean water in the tank and drain through the used tip by switching the appliance on.

- Appliance disinfection. Pour 200-300 ml of 40 % medical alcohol solution into the tank and drain it through the tips by switching the appliance on. One hour later, rinse the tank with clean water, pour the filtered clean water in the tank and drain it through the tips.

STORAGE Store the assembled appliance in cool and dry place, far from children.

TROUBLE SHOOTING

Trouble	Purpose	Trouble shooting
Leakage between the tip and rotary mechanism of the tip.	The tip was installed incorrectly.	Detach the tip and reinstall it.
Tooth pocket tip. Leakage of the soft end.	Soft end is torn.	Replace the tip.
Hose leakage.	Hose is damaged.	Replace hose in the service center.
Water tank leakage.	Water tank installed on the base incorrectly.	Reinstall the water tank.
Low water pressure in the appliance / water does not flow.	The water tank is installed onto the base incorrectly. No water in the tank.	Reinstall the water tank, then fix it in the base by slightly pressing. Fill the water tank.
The appliance does not turn on.	The motor does not work.	Check the presence of current in mains.

PRECAUTIONS

• the appliance is for home use only. Do not use for industrial

purposes. Use the appliance strictly in accordance with your dentist recommendations. • before the first operating make sure that the voltage in your wall outlet is the same as on the appliance's rating label marked. • if you dropped the appliance and it has some visible defects, discontinue using it, bring to an authorized service center for checking. • do not immerse the appliance or power cord into water. • do not touch the power cord with wet hands. • instruct the children or people with physical capabilities to use the appliance correctly. • do not use the appliance on a damp surface/floor. • do not let working appliance be unattended. • always unplug the appliance before cleaning. • do not use any abrasive, chemical or alcohol-contained cleaners, polish for cleaning the appliance. • do not use the damaged appliance or appliance with damaged cord/plug. • do not let the cord hang over the edge of table or touch any hot surfaces. • use only enclosed tips. • unplug the appliance when it is not in use. • use appliance only for its direct purpose. • do not expose the appliance to fall down or to be kicked. Repair the appliance only in an authorized service center. • Do not repair the appliance by yourself! If the opening of the appliance is detected, the appliance's guarantee will be canceled. • while using the appliance, please, observe the safety precautions. • incorrect operation and improper handling can lead to malfunction of the appliance and injuries to the user. • do not use the appliance while taking bath/shower. • do not immerse appliance in water. • if the appliance fell down into the water, unplug it immediately and do not use it. Bring it to an authorized service center for checking. • do not direct water jet under your tongue, into the ear, nose or at any body parts. Water jet pressure can cause serious harm to these body parts. • do not fill water tank with water, if its temperature exceeds 40°C. This may cause burns of mucous coats. • fill the water tank with water or fluids (rinsing or therapeutic fluids), recommended by your dentist and nothing else. The liquid should not leave sediment or be viscous, this may soil the appliance. • in case you are having oral treatment, or suffer with severe paradontosis consult the dentist before using it. • do not insert foreign objects into the appliance. • do not let power cord and appliance touch hot surfaces. • do not use bleaching or iodine-containing solvents. • it is necessary to pour water remnants out of water tank after use. • do not use this appliance if you have

jewelries in the mouth; take them away before using the appliance.

- do not use the appliance if you have open sore on tongue/ in the mouth.
- give a high attention when the appliance is being used near the children. Do not let children use or play with the appliance.

given appliance is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience or knowledge, unless they have been given supervision or instruction concerning use of appliance by a person responsible for their safety.

- do not use this appliance for more than 5 minutes in each two-hour period. Longer usage can overheat the motor and cause appliance damage as well as void the warranty.

USING MEMO Attention! Fill the water tank with filtered water only! Usual tap water contains salt, sand and lime admixtures that can quickly soil the appliance. The appliance that comes to the service center with clog caused by dirty water, bad appliance cleaning; or using soiling liquids, the appliance's guarantee will be canceled for incorrect handling.

- After using some mouth rinse in this appliance, empty the water tank and clean the appliance. Rinse water tank with warm water and fill with filtered water. Run water through the used tip to clean the water supply hose and the used tip. For this plug the appliance in, choose any water pressure level, switch the appliance on and perform the appliance cleaning over the sink.

Tips of appliance do not have guarantee, as they are considered to be personal hygiene items.

DISPOSAL Dispose the appliance in an environmentally friendly manner. It should not be disposed with the normal household waste. Please dispose it at a recycling centre for used electrical and electronic devices. You can obtain more information from your local authorities.

NOTE: The manufacturer reserves the right, without prior notice, to make changes in the product that do not impair its technical characteristics, and are the result of improvements in its design or production technology.

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