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## DENTAL WATER FLOSSER

MODEL: VIP-009-B

Dear customer! Congratulations with the purchase of home appliances under VES Electric brand. From our side we guarantee excellent quality, perfect functioning and safety of your appliance upon keeping to all the rules of home exploitation. Before using the appliance read the instruction manual carefully.

### Technical data:

Voltage: 100-240 V Frequency: 50/60 Hz Power: 18 W

Maximum pressure of the jet: 879 kPa

Capacity of water tank: 600 ml

### Recommendations for buyer

When buying the appliance, check it for the absence of mechanical damage of the construction, check the external design and packing, the functioning, the presence of all accessories in set. Note: operation checking - is the duty of the seller. Please check the warranty card and the accuracy of its filling (presence of the seller's seal). You must save the bill, instruction manual, warranty card until the end of the warranty period. In case of bill and warranty card loss they cannot be renewed. In the case of contacting the service center, if it is impossible to determine the date of product purchasing, the warranty period is determined from the manufacturing date.

### DEVICE DESCRIPTION

1. Water tank cover (there is a compartment for tips storage)
2. Periodontal nozzle
3. Water tank
4. The mechanism of tip rotation
5. Tip eject button
6. Handle
7. Pause button

8. 0/1 button
9. Pressure regulator
10. Handle place
11. Body
12. Suction feet
13. Tips jack


### BEFORE FIRST USE

Remove the appliance from the package. Do not leave the packing (plastic bags, polystyrene etc.) within the

---

reach of children since they could be dangerous. Keep the package for appliance storage. After transportation in low temperatures leave the appliance in the room for 2 hours before use.

#### Assembly scheme

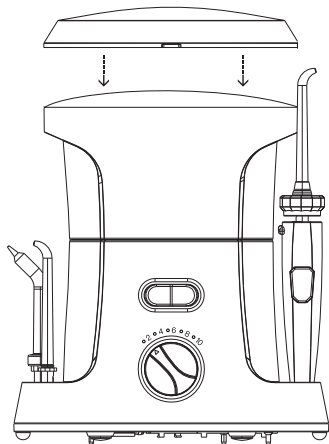
#### PREPARATION FOR USE

Detach the water tank from the base.

**There is a rubber plug at the bottom of the water tank. Before the operation make sure the rubber plug is in its place and inserted correctly.**

Fill the water tank with warm water. Don't fill water tank with water, if its temperature exceeds 40°C. This may cause burns of mucous coats. Place the water tank back to the base, slightly press it for tight fitting and close with the lid.

**Note:** If your physician recommended you to take antibiotics

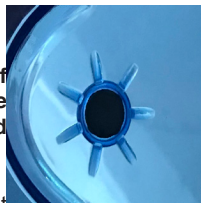---

before having oral hygiene procedures, consult your dentist before using given appliance or any other appliance for oral hygiene.

#### Water tank filling

**Use only the filtered water! Tap water contains salt, sand, lime, etc., which will clog the appliance. The clogged / dirty unit or the unit with poor cleaning voids the warrantee.**

Fill water tank with warm water, close the lid.

**Attention!** Do not use bleaching or iodine-containing fluids.

Fill the water tank with water or fluids (rinsing or medicinal fluids), recommended by your dentist and nothing else. The liquid should not leave sediment or be viscous, this may soil the appliance.

If you used fluid, you should pour out its remnants and clean the appliance. For this fill the water tank with filtered water and drain water through used tip, turning

#### Attention!

Contact your dentist if:

- Do you doubt the harmlessness of the use of the device.
- After using the device appears excessive bleeding of your gums, or if bleeding continues after 2 weeks of using the device.

#### Tip placement

Place the tip in the handle until the click.

#### Tip detaching

To detach the tip press on the tip eject button and take the tip out.

**Attention!** Do not press the tip eject button / change the tip while the device is working.

Do not replace the tips too often, as the holes and seal mechanism are subject to wear.

**Tips of the device, as the subject of personal hygiene, do not fall under the guarantee.**

## Water pressure regulation

Set the water pressure regulation to the left side is for minimum water pressure and gradually set it right for higher pressure.

## Location in the mouth


Locate the tip in the mouth. Direct the tip end at the teeth. Switch the appliance on. Set the required water level. Direct water jet, turning the appliance. Close your lips to avoid water splashing, but do not close your mouth allowing water freely flow away. For better result please start with side teeth and molars, moving it to fore teeth. Move tip from tooth to tooth making pauses between moves. Go on cleaning until teeth surrounding zones and zones between teeth are clean.

## Using tips


**Caution!** Do not direct water jet under your tongue, into the ear, nose and any body parts except mouth. Water jet pressure can cause serious harm to these body parts. Do not use this appliance if you have jewelleries in the mouth; take them away before using the appliance. Do not use the appliance if you have open sore on tongue / in the mouth.

**Pause** You can temporarily stop the irrigation by pressing the Pause button on the handle.

## TIPS

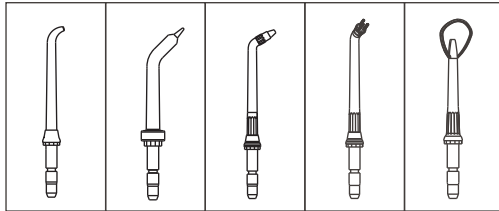Jet tip      Periodontal tooth pocket tip      Orthodontic tip      Plaque seeker tip      Tongue cleaner
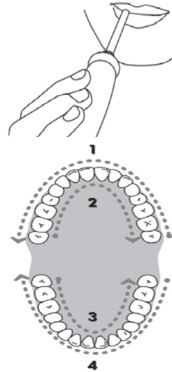
## Standart jet tip

Jet tip is used for regular teeth cleaning.

## Tongue cleaner

It is used for cleaning the tongue from plaque and bacteria. Place the tip on the center of tongue. Use some slight force to pull the tip to the direction of the lips. Repeat several times to clean the entire surface of the tongue. In the process it is possible to gradually increase the pressure of the water if it does not bring discomfort.


## Plaque seeker tip

It is used for cleaning crowns, implants, bridges. The effect is achieved due to the unusual shape of the trident-shaped nozzle.


## Periodontal tooth pocket tip

Tip allows supplying water or antibacterial solvent directly into the teeth pockets. Direct soft end of the tip at the tooth at 45 degrees angle and gently insert the tip under the gum line, into the tooth pocket. Turn the appliance on and perform cleaning moving the tip between the teeth and along the gum line.

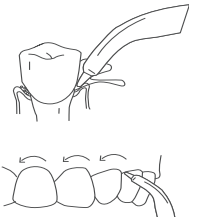**Orthodontic tip** is used for cleaning inaccessible places: around dental brackets, crown works, frameworks and dentures. Direct the tip at the teeth at 90 degrees angle, making slight effort move the tip from molars to the fore teeth clean the space along gum line, hold the tip between teeth for some time to clean the space between them and around dental bracket, then clean inner side of the tooth and get to the next tooth.
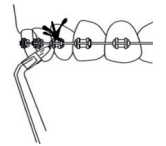
## CLEANING

Unplug the appliance from mains. Detach water tank. Pour water remnants out of water tank. Rinse water tank with clean water, let dry. Wipe the base with wet cloth. Install water tank back onto the base. Roll power cord carefully. For appliance cleaning use only the clean water or neutral cleansers. Do not use the caustic or abrasive cleansers (e.g. vinegar or decalcifying agents) as they may damage the appliance. Do not wash the tips with water hotter than 40°C.

- **If you used the mouth wash** its remnants must be cleaned away. For this clean the tank with clean water, pour the filtered clean water in the tank and drain through the used tip by switching the appliance on.

- **Appliance disinfection.** Pour 200-300 ml of 40 % medical alcohol solution into the tank and drain it through the tips by switching the appliance on. One hour later, rinse the tank with clean water, pour the filtered clean water in the tank and drain it through the tips.

## Storage

Keep the unit assembled in a clean, dry place, away from children.

When transporting and storing the device, observe the following rules:

- The unit must be empty and completely dry.
- After transportation of device at low temperatures must be left in the room for at least 2 hours before using.

## Problems and solutions

Problem	Possible reason	Solution
Leakage between the tip and the tip rotating mechanism.	Tip incorrect installation.	Eject the tip and install correctly.

Leakage of periodontal tooth pocket tip end.	Tip if torn.	Change the tip.
Leakage of hose.	Hose is torn.	Change the hose at service center.
Leakage of water tank.	Incorrect installation of the water tank or there is no rubber plug.  Rubber plug placed incorrectly.	Install the water tank correctly. Install the rubber plug back. If plug lost, address to service center. Install the rubber plug correctly.
Low water pressure or there is no water coming.	Water tank is incorrectly placed on the base. There is no water in tank.	Place water tank well.  Fill tank with water.
The appliance is switched on.	Motor is not working.	Make sure there is current in net.
Tips and tank are correctly installed. There is water in tank but the water is not running.	There may be the airlock in the device.	To eliminate the airplug turn the unit on and gently retract the air from a tip without swallowing suddenly entered the water.

## PRECAUTIONS

**Use only the filtered water! Tap water contains salt, sand, lime, etc., which will clog the appliance. The**

---

**clogged / dirty unit or the unit with poor cleaning voids the warrantee.**

- the appliance is for home use only. Do not use for industrial purposes.
- use the appliance strictly in accordance with your dentist recommendations.
- before the first operating make sure that the voltage in your wall outlet is the same as on the appliance's rating label marked.
- if you dropped the appliance and it has some visible defects, discontinue using it, bring to an authorized service center for checking.
- do not drop or kick the appliance.
- do not immerse the appliance, power cord or plug in water.
- people that can not operate the appliance themselves, or people with low oral sensability should not use this appliance.
- keep the appliance far from children, breaking of this rule may lead to an accident and dangerous situations like tips swallowing.
- instruct and look after the children or people with physical capabilities to use the appliance correctly.
- do not let working appliance be unattended.
- do not use any abrasive, chemical or alcohol-contained cleaners, polish for cleaning the appliance.
- do not use the damaged appliance or with damaged cord/ plug.
- do not let the cord hang over the edge of table or touch any hot surfaces.
- use only enclosed tips.
- use the appliance for the direct purpose only.
- repair the appliance only in an authorized service center. Do not repair the appliance by yourself! If the

---

opening of the appliance is detected, the appliance's guarantee will be canceled.

- while using the appliance, please, observe the safety precautions.
- incorrect operation and improper handing can lead to malfunction of the appliance and injuries to the user.
- do not leave the appliance where it can drop in water.
- if the appliance fell in the water, discontinue it's use  
Bring it to an authorized service center for checking.
- do not direct water jet under your tongue, into the ear, nose or at any body parts. Water jet pressure can cause serious harm to these body parts.
- do not fill water tank with water, if its temperature exceeds 40°C. This may cause burns of mucous coats.
- fill the water tank with water or fluids (rinsing or medicinal fluids), recommended by your dentist and nothing else. The liquid should not leave sediment or be viscous, this may soil the appliance.
- in case you are having oral treatment, or suffer with severe paradontosis consult the dentist before using it.
- do not insert foreign objects into the appliance.
- do not use bleaching or iodine-containing solvents.
- it is necessary to pour water remnants out of water tank, then leave the it to fully dry after use otherwise mold and bacteria can badly affect your health, lead to unpleasant smell of the appliance and worsen appliance's work.
- do not turn on appliance with empty water tank.
- don't use the appliance in the vicinity where explosive gas, paint thinners, gasoline, paint or other flammable liquids are used or stored.
- do not use this appliance if you have jewelries in the mouth; take them away before using the appliance.
- do not use the appliance if you have open sore on

---

tongue / in the mouth.

- give a high attention when the appliance is being used near the children. Do not let children use or play with the appliance.
- do not use the appliance outdoors.
- given appliance is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience or knowledge, unless they have been given supervision or instruction concerning use of appliance by a person responsible for their safety.

**Tips of appliance do not have guarantee, as they are considered to be personal hygiene items.**

**DISPOSAL** Dispose the appliance in an environmentally friendly manner. It should not be disposed with the normal household waste.

■ Please dispose it at a recycling centre for used electrical and electronic devices. You can obtain more information from your local authorities.

**NOTE: The manufacturer reserves the right, without prior notice, to make changes in the product that do not impair its technical characteristics, and are the result of improvements in its design or production technology.**

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